
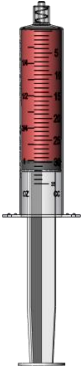
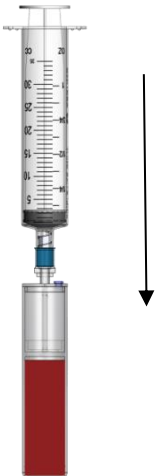
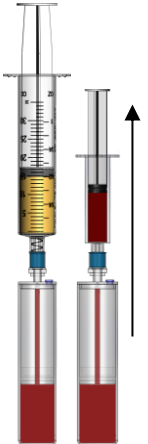


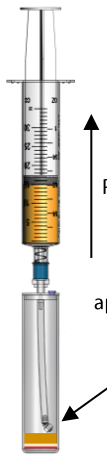


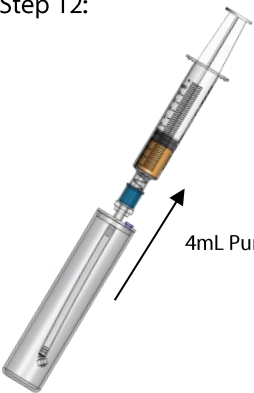


**GS30-PURE II PROTOCOL B \*\* PLEASE DISCARD RED VENTED CAP FROM CONCENTRATING DEVICE BEFORE USE \*\***

Note: Always swab self-sealing port with sterile alcohol prior to accessing with a sterile syringe

<p>Step 1:</p>  <p>Draw 3mL of Citrate Anticoagulant into 30mL Syringe</p>	<p>Step 2:</p>  <p>Draw 25mL whole blood from the patient, filling the syringe to 30mL</p>	<p>Step 3:</p>  <p>Load anticoagulated whole blood into the <b>Concentrating Device</b></p>	<p>Step 4:</p> <p>Counterbalance and process the <b>Concentrating Device</b> at</p> <p><b>Executive Series Centrifuge</b> 1 minute 4200 RPM</p> <p><b>Euro Series Centrifuge</b> 2 minutes 3500 RPM</p>
<p>Step 5:</p>  <p>Using the 30mL syringe, aspirate the platelet plasma suspension (PPS) until RBC fills the aspirating pipe.</p> <p>Then using the 3mL syringe aspirate 1mL of RBC</p>	<p>Step 6:</p>  <p>Transfer the platelet plasma suspension (PPS) and the RBC into the <b>Concentrating Accessory</b></p>	<p>Step 7:</p> <p>Counterbalance and process the <b>Concentrating Accessory</b> at</p> <p><b>Executive Series Centrifuge</b> 5 minutes 4200 RPM</p> <p><b>Euro Series Centrifuge</b> 6 minutes 3500 RPM</p>	<p>Step 8:</p>  <p>Platelet Concentrate Buffycoat</p> <p>Platelet concentrate and neutrophil cytokine buffycoat separates out at the bottom of the <b>Concentrating Accessory</b></p>
<p>Step 9:</p>  <p>Aspirate Platelet Poor Plasma</p> <p>Leave approximately 4mL</p> <p>Aspirate platelet poor plasma from the <b>Concentrating Accessory</b> Leave 4mL of plasma.</p>	<p>Step 10:</p>  <p>Attach the 12mL syringe and swirl to resuspend the platelet buffycoat into the plasma.</p>	<p>Step 11:</p>  <p>Tilt to immerse Aspirating Pipe into PurePRP®</p> <p>Tilt to immerse the Aspirating Pipe into the PurePRP®</p>	<p>Step 12:</p>  <p>4mL PurePRP®</p> <p>Extract the PurePRP® into the 12mL syringe.</p>